

GENESIS ALEX

“The most complete and technologically
Advanced solution in SKIN TONING & SKIN LIFTING.”

Alexandrite laser
advantages
+
Without
consumables



Integrating the most advanced technology in dealing with skin lifting and pigmentation, GENESIS ALEX guarantees the highest level of effectiveness & safety during treatment.

This incomparable benefit makes GENESIS ALEX as the only alexandrite laser with zero treatment cost, giving the doctor the advantage in challenging with the ongoing competition in pigmentation and lifting fields.

At the same time, GENESIS needs practically no consumables, which allows the doctor to ensure zero based treatment cost.

GENESIS ALEX

WAVELENGTH	755nm
PULSE DURATION	5ms ~70ms (Adjustable)
CONTROL	10.4" touch screen LCD
CONTAC TIP	Sapphire contact chill tip
SPOT SIZE	15mm x 15mm
MAX REPETITION RATE	1Hz ~ 10Hz
SKIN COOLING	Semiconductor method (0°C ~5°, adjustable)
DIMENTION	350(W)x400(L)x950(H) mm
WEIGHT	50kg



2F, 325-2, Ttukseom-ro, Seongdong-gu,
Seoul, Korea 04779
T: +82-2-465-2990~1
+82-2-499-5573
F: +82-2-465-2911
ilaserlab@ilaserlab.com
ilaserlab@naver.com
www.ilaserlab.com



ALEXANDRITE
REBORN

GENESIS ALEX




The absolute GAME CHANGER FOR SKIN LIFTING & SKIN TONING


GENESIS ALEX is a new, technologically advanced ALEX Laser for skin lifting & pigmentation, which combines conventional alexandrite laser with the advanced technology of diode laser.

The wavelength of 755nm in GENESIS Alex enhances the energy absorption rate in the melanin, and its large spot with integrated cooling system, guarantees even faster, painless and effective treatment.

GENESIS ALEX is offering a complete and technologically advanced solutions in treatment of superficial pigmentation, ensuring the highest level of effectiveness and safety for every skin photo-type.

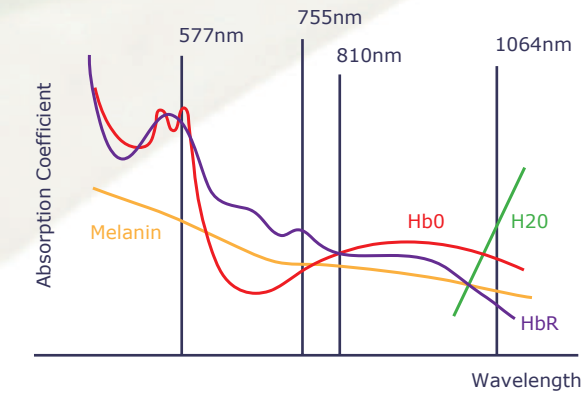

ALEX 755nm wavelength:
Ideal wavelength for pigmentation

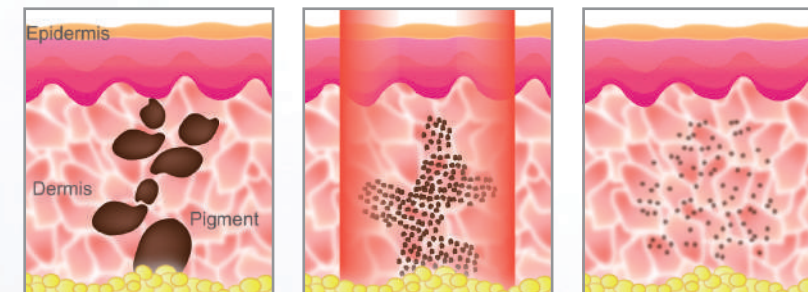

GENESIS ALEX laser system can successfully treat a wide range of cosmetically undesirable benign pigmented lesions without damaging adjacent tissues around melanin, by emitting of trains of micro sec unit ON-TIME 755nm wavelength

GENESIS ALEX integrates the most advanced and patented technology and guarantees the highest level of safety and performance during treatment

SPEED UP TO 10HZ
Saves both doctor's and patient's valuable time.

ANY KIND OF PHOTOTYPE, ANY TIME IN YEAR

Safe treatments for all phototypes, all year round, due to gradual heat accumulation technology.

PRACTICALLY PAINLESS TREATMENT
In-motion technology combined with Alexandrite.

CONTACT COOLING
maintains skin temperature at desired levels throughout the whole treatment.

TOTAL COVERAGE
With large spot size of 15mm x 15mm, GENESIS ALEX covers every targeted area during treatment.




GENESIS ALEX

SYSTEM ADVANTAGES


TRAINS OF PULSE IN PULSE WAVE

GENESIS Alex emits trains of PIPW to give an effective thermal accumulation in follicles without damages in skin tissue.

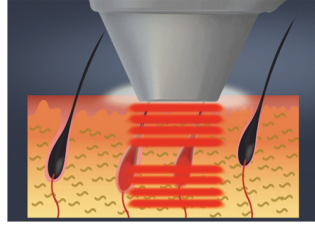


DEPTH OF PENETRATION

GENESIS Alex realizes largest square 15mm x 15mm spot size. Larger spot size enables deeper penetration, and it leads to improved absorption rate of energy at the target follicle.

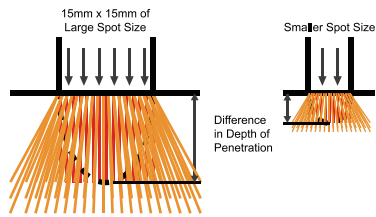

BENIGN PIGMENTED LESIONS

GENESIS Alex can successfully treat a wide range of cosmetically undesirable benign pigmented lesions without damaging to adjacent tissues around melanin by emitting of trains of micro sec unit pulses.

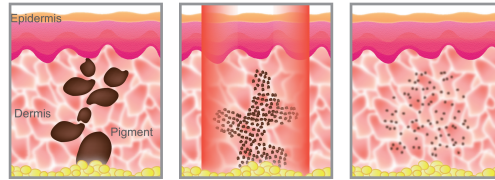

***Semi Permanent Module
(2 Year warranty)**

***Low maintenance cost for ownership**

Convenient **GUI SYSTEM**
GUI system for easy use

